


# IMPLANT ORAL SURGERY HANDPIECES

## Low Speed Handpieces for Implantology & Oral Surgery

*For implantology, ITL Dental  
offers handpieces for surgical  
drilling and bone tapping*

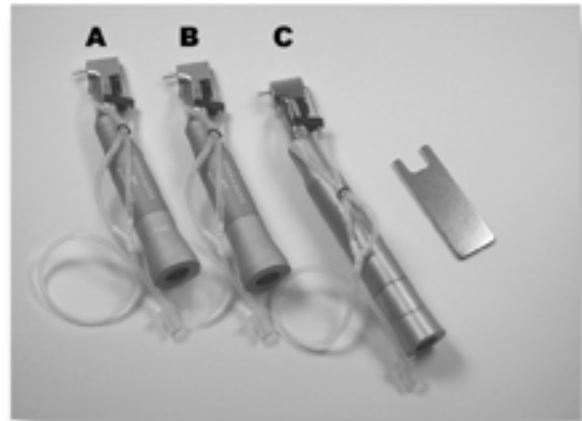

- A** The 18:1 Reduction Latch Contra Angle handpiece for surgical drilling is the only handpiece with a sealed latch head that requires no lubrication or cleaning maintenance. Simplify cleaning procedures for your staff and reduce costly handpiece maintenance requirements. The handpiece has internal and/or external irrigation. Product Code: 28181-S
- B** The 32:1 Reduction Latch Contra Angle handpiece has internal and/or external irrigation. This handpiece can be used for both surgical drilling and bone tapping for selected surgical drilling systems. Product Code: MD5063
- C** The 400:1 Reduction Latch Contra Angle comes with a heavy-duty latch head to absorb the high torques of bone tapping and mechanical implant insertion. Internal and/or external irrigation can be provided from the latch head. Product Code: 28402



## Oral Surgery Hand pieces: High Speed or High Torque

- A** The 1:3, 90,000 rpm straight surgical handpiece provides external irrigation during use. It is available with the "E" type connection- (Product Code 28213) for electric systems, air motor converters or 2-hole air connection. An extra bur guard is included with the handpiece. Product Code 28213 ("E-type") or Product Code 28113 (with air connection.)
- B** The 1:1, 30,000 rpm straight surgical handpiece provides external irrigation during use. It has the "E" type connection for use with electrical systems or air motor converters. An extra bur guard is included with the handpiece. Product Code: 28211
- C** Replacement surgical handpiece bur guards. Product Code: 28552

**ITL DENTAL**

31 Peters Canyon  
Irvine, California 92606  
PH: 949.223.8950 FX: 949.223.8960  
**800.277.0073**  
www.itldental.com sales@itldental.com

*"Innovative Surgical Handpieces"*